

Lunar® U PRO

Intelligent Powered Stapler

The innovative stapling technology that combines compression and Varied-Height-Staple technology.

Lighter, Power, Smarter

Real-time detection of tissue thickness and automatic stepless adjustment, to improve surgical and patient outcomes.

FEATURES


Comfort and Control

Lunar® U Pro delivers the ergonomic comfort and control you need to perform at your best. The OLED screen displays real-time feedback during the stapling process creating more precise mechanical visualisation and control. Bidirectional articulation provides access to difficult areas, to improve surgical visualisation and control. The detachable “Sled” can assist in the accessing and efficient handling of blood vessels.

- ✓ Smarter Powered Technology
- ✓ Stepwise Varied-Height-Staple
- ✓ Predictable & Reproducible Result
- ✓ Positive Patient Outcomes

Lunar® U PRO Powered Stapler

RE-ORDER	JAW LENGTH	WORKING LENGTH	CARTON QUANTITY
DSAS45	45mm	280mm	3 pcs / Carton
DSAS60	60mm	280mm	3 pcs / Carton
DSAM45	45mm	340mm	3 pcs / Carton
DSAM60	60mm	340mm	3 pcs / Carton
DSAL45	45mm	440mm	3 pcs / Carton
DSAL60	60mm	440mm	3 pcs / Carton

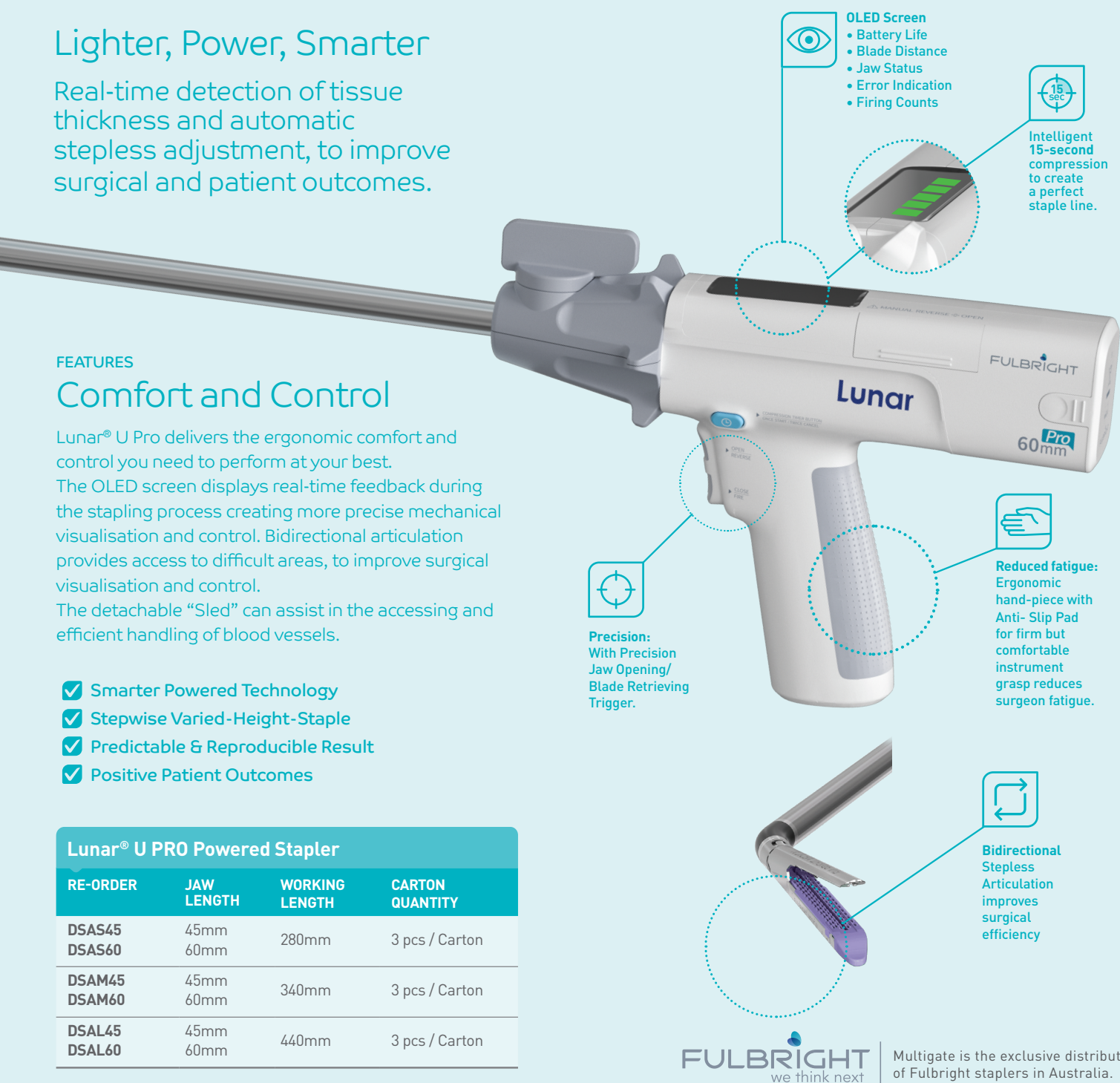

Lunar® U PRO

Stepwise Staple

Varied-Height-Staple technology.

Maximised Sharpness Guaranteed

Intelligent functionality with a new blade built into each reload, ensuring sharpness every time.

FEATURES

Improved Patient Outcomes

3 rows of staples (each of varying heights) reduce tissue stress and accommodates different tissue structure.

Improved Perfusion Volume



Varied-Height Staples deliver high perfusion volume and might deliver better wound healing and clinical outcomes.

Lunar® U PRO Powered Stapler

RE-ORDER	JAW LENGTH	COLOUR	ROW	OPEN STAPLE HEIGHT (mm)	CARTON QUANTITY
DACW45	45mm	White	6	2.5	12pcs / Carton
DACW60	60mm				
DACT45	45mm	Tan	6	2.5 / 3.0 / 3.5	12pcs / Carton
DACT60	60mm				
DACP45	45mm	Purple	6	3.5 / 3.8 / 4.1	12pcs / Carton
DACP60	60mm				
DACX45	45mm	Black	6	4.1 / 4.4 / 4.7	12pcs / Carton
DACX60	60mm				


B2B Staircase Deck

- Back-to-back *3 varied staple staircase structure in height.
- Accommodates tissue structure
- Prevents unexpected bleeding.

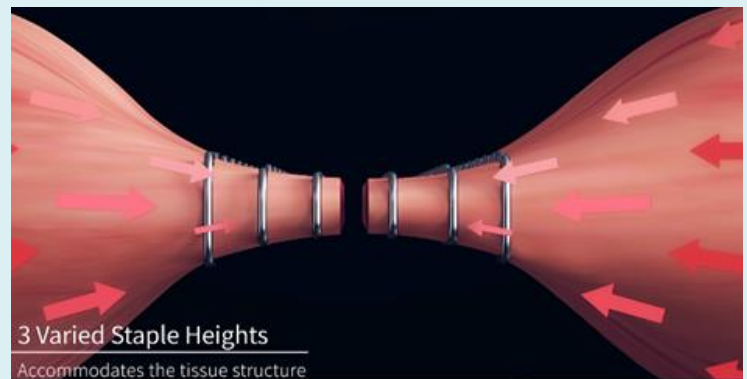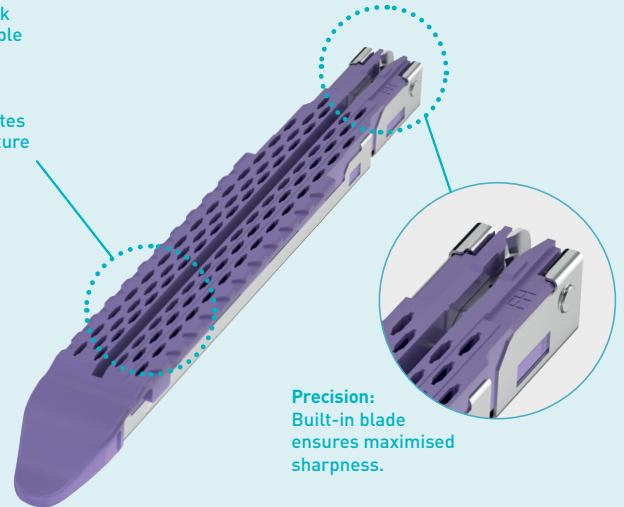

Optional Accessory

Curved Tip: Improving the efficiency of handling blood vessels.

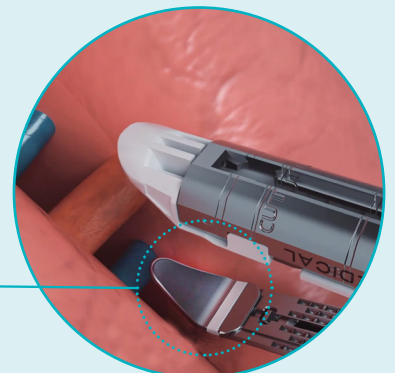

FULBRIGHT
we think next

Multigate is the exclusive distributor of Fulbright staplers in Australia.

MULTIGATE

For more information or to book a trial, contact Multigate Territory Manager: (02) 8718 2888 or visit multigate.com.au

MG Lunar® U Pro v4.0 14JAN2026