

MULTIGATE

Pro-Zone[®]



High Performance Sterilisation Wrap
Innovation Focus: Improved Barrier Performance

V 2.0



At Multigate we combine our deep understanding of how you work with our design capabilities to innovate and better serve you, your patients and our partners.

For 30 years we've created medical products, services and experiences that give you the advantage to perform at your best. We continue to do this.

Inspiring performances and better patient outcomes that will lead the world of hospital care to new standards.

We are Multigate.

Addressing CSSD clinician concerns about the time, costs and effectiveness of sterilising instruments, we introduced our simple to use, cost saving Pro-Zone® range of single use sterilisation wraps. Made from polypropylene and using our 5 layer SSMMS protection technology, Pro-Zone® wraps offer an economic and time efficient sterilisation solution.



INFECTION PREVENTION

Pro-Zone®’s strong 5 layer construction reduces the risk of tears, punctures and abrasions, and delivers a low linting, high performing sterile barrier system.

- ✓ Infection prevention
- ✓ Process efficiencies
- ✓ Comfortable handling
- ✓ Space saving
- ✓ Recyclable
- ✓ Low linting
- ✓ Anti-static

PERFORMANCE

Streamlined Process.

Our Pro-Zone® range eliminates the costly decontamination processes associated with reusable sterilisation systems and delivers faster instrument processing.

The fabric is soft, has low memory and provides excellent draping capability allowing Pro-Zone® wrap to double as a sterile field, permitting the aseptic removal of pack contents and further streamlining theatre processes.

In-Service Support.

We understand that every CSSD is different. Our Sterilisation Specialists are trained in CSSD principles and practises. They are available to:

- Provide advice
- Deliver in-service training and support
- Evaluate sterilisation wrap performance.

FEATURES

Comfortable Handling.

Pro-Zone® wraps are anti-static, easy to fold, gentle on your hands and light offering a packaging system that is easy to handle.

Space Saving.

Pro-Zone®’s compatibility with a broad range of sterilants and its ability to cater for all instrument shapes and sizes with just a few different weights and specifications, reduces the number of sterilisation components and returns valuable space to CSSDs.

Recyclable.

The entire Pro-Zone® range is recyclable.



Microbial barrier that delivers improved barrier performance



Easy to fold material with low memory to provide excellent draping capability



TECHNOLOGY



Our Pro-Zone® range is made from polypropylene, utilizing state of the art 5 layer (Spunbond-Spunbond-Meltblown-Meltblown-Spunbond) technology.

DESIGN

Our entire Pro-Zone® range of sterilisation wraps are a versatile high quality wrapping system that delivers many performance advantages including:

- Strength and resistance to punctures, tears and abrasions*
- Anti-static coating reduces static “shocks” to clinicians
- Low linting
- Soft and easy to fold
- Excellent drapability and low memory
- Excellent microbial barrier
- Effective sterilant penetration

* when correct handling techniques are used

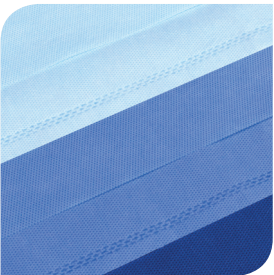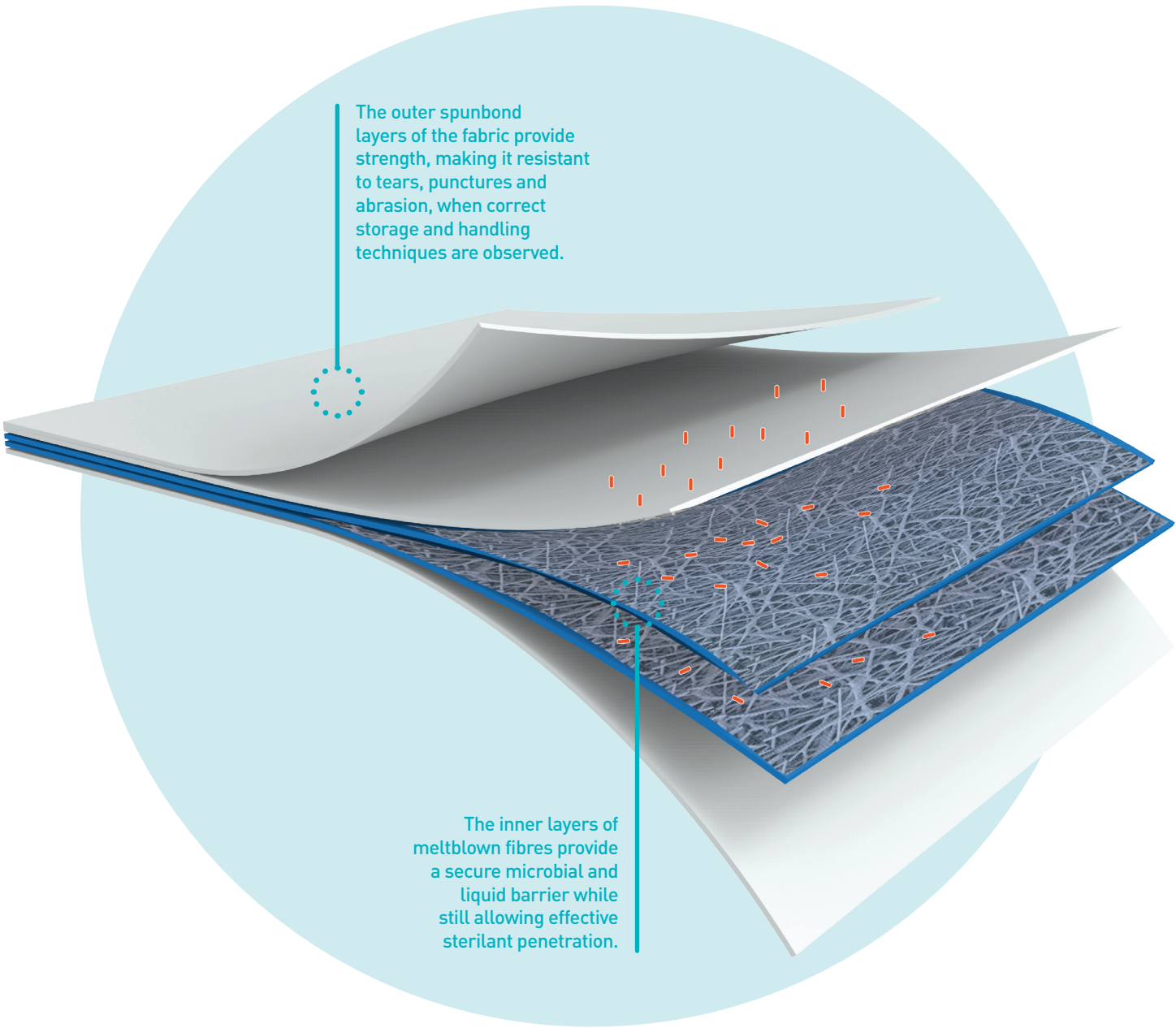

Pro-Zone®
The original Pro-Zone® sterilisation wrap.

- Weights:**
- 47gsm
 - 61gsm
 - 75 gsm

- Compatible with:**
- Steam
 - ETO
 - Low temperature gas plasma
 - Low temperature steam/formaldehyde

- Standards Compliance:**
- AS/NZS 4187:2014
 - ISO 11607



Pro-Zone®+2
Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.

- Weights:**
- 47gsm
 - 61gsm
 - 75 gsm

- Compatible with:**
- Steam
 - ETO
 - Low temperature gas plasma
 - Low temperature steam/formaldehyde

- Standards Compliance:**
- AS/NZS 4187:2014
 - ISO 11607



Pro-Zone®+4G
Pro-Zone®+4G offers the protection of 2 bonded sheets, plus the added advantage of an absorbent layer to wick any residual condensate that may compromise the integrity of the wrap.

- Weights:**
- 61gsm
 - 75 gsm

- Compatible with:**
- Steam
 - ETO
 - Low temperature steam/formaldehyde

- Standards Compliance:**
- AS/NZS 4187:2014
 - ISO 11607



Pro-Zone®

The original Pro-Zone® sterilisation wrap.



Pro-Zone®+2

Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.

Pro-Zone® (47) COLOUR: BABY BLUE			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-465	30 x 30	47	1000
30-464	45 x 45	47	500
30-460	61 x 61	47	500
30-470	71 x 71	47	300
30-461	91.5 x 91.5	47	150
30-467	100 x 106	47	150
30-466	100 x 140	47	100
30-462	121.5 x 121.5	47	100
30-468	130 x 152.5	47	100
30-463	183 x 137	47	100


Pro-Zone® (61) COLOUR: BLUE			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-471	71 x 71	61	200
30-472	91.5 x 91.5	61	150
30-475	106 x 120	61	50
30-476	106 x 140	61	50
30-473	121.5 x 121.5	61	50
30-477	130 x 152.5	61	50
30-478	137 x 137	61	50
30-474	137 x 183	61	50

Pro-Zone® (75) COLOUR: DARK BLUE			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-484	106 x 120	75	50
30-483	106 x 140	75	50
30-480	121.5 x 121.5	75	50
30-481	137 x 183	75	40

Pro-Zone®+2 (47) COLOUR: BABY BLUE / LILAC			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-551Z	30 x 30	47	480
30-552Z	45 x 45	47	480
30-553Z	60 x 60	47	240
30-554Z	76 x 76	47	144
30-555Z	91 x 91	47	144
30-556Z	101 x 101	47	120
30-557Z	121 x 121	47	48
30-558Z	137 x 137	47	48
30-559Z	137 x 182	47	48

Pro-Zone®+2 (61) COLOUR: BLUE / LILAC			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-560Z	71 x 71	61	96
30-561Z	91 x 91	61	72
30-562Z	101 x 101	61	72
30-563Z	121 x 121	61	24
30-564Z	137 x 137	61	24
30-565Z	137 x 182	61	24

Pro-Zone®+2 (75) COLOUR: DARK BLUE / LILAC			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-570Z	91 x 91	75	72
30-571Z	101 x 101	75	48
30-572Z	121 x 121	75	24
30-574Z	137 x 137	75	24
30-573Z	137 x 182	75	24

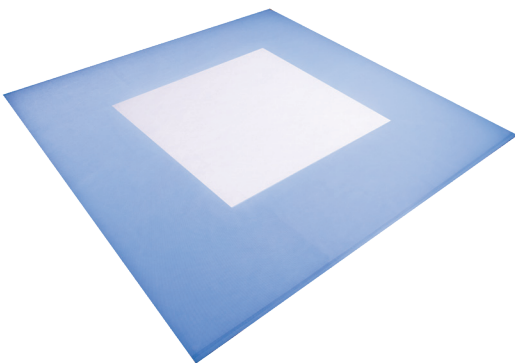

Pro-Zone®+4G

Pro-Zone®+4G offers the protection of 2 bonded sheets, plus the added advantage of an absorbent layer to wick any residual condensate.

Pro-Zone®+4G (61) COLOUR: BLUE			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-863	121 x 121cm (absorbent layer 70 x 70)	61	24
30-864	130 x 152.5cm (absorbent layer 70 x 70)	61	24

Pro-Zone®+4G (75) COLOUR: DARK BLUE			
RE-ORDER	SIZE (CM)	WEIGHT (GSM)	CARTON QTY (PIECES)
30-872	121 x 121 (absorbent layer 70 x 70)	75	24
30-874	137 x 137cm (absorbent layer 70 x 70)	75	24
30-873	137 x 182cm (absorbent layer 70 x 70)	75	24



Towels and Tray Liners

Made from non-woven spun lace cellulose/ polyester fabric, MediClean™ is ideal for a wide range of dry or wet CSSD applications. MediClean™ towels are ideal for use as a tray liner and applications where the towel is soaked in water for long periods of time.

MediClean				
RE-ORDER	SIZE (CM)	COLOUR	PIECES/BOX	BOXES/CARTON
57-112	31.8 x 33 (small)	○ White	100	12
57-113	60 x 35 (large)	○ White	100	12



Protective Corner

Our protective corner helps to reduce piercing risk while wrapping.

Protective Corner		
RE-ORDER	COLOUR	CARTON QTY (PIECES)
DP 2704	● Orange	300

Empowering
performance
since 1986

Multigate Medical Products

27 Llewellyn Avenue

Villawood NSW 2163

Australia

T +61 (0)2 8718 2888

F +61 (0)2 8718 2889

enquiries@multigate.com.au

www.multigate.com.au

For more information please contact your
local Business Development Manager:

(02) 8718 2888 or visit **multigate.com.au**