

# MULTIGATE SafeCount™



## Needle Counter 20 count / 40 capacity

STERILE

PRECISION,  
ORGANISATION  
AND PROTECTION


The SafeCount™ needle counter securely manages and counts sharps with innovative features designed for precision, organisation, and enhanced protection.

RE-ORDER CODE

CEA-0010S

### FEATURES AND BENEFITS

- ✓ Notches created in unit for instrument organisation
- ✓ Novel foam slots assist with suture management
- ✓ Easy blade lift from undulating magnetic strip
- ✓ Printed numbers to assist with counting sharps
- ✓ X-ray detectable secure lock ensures needle counter is closed and ready for safe disposal
- ✓ Adhesive fixation on rear to keep needle counter in place
- ✓ Easy to identify product has been coloured yellow, the industry standard colour for sharps

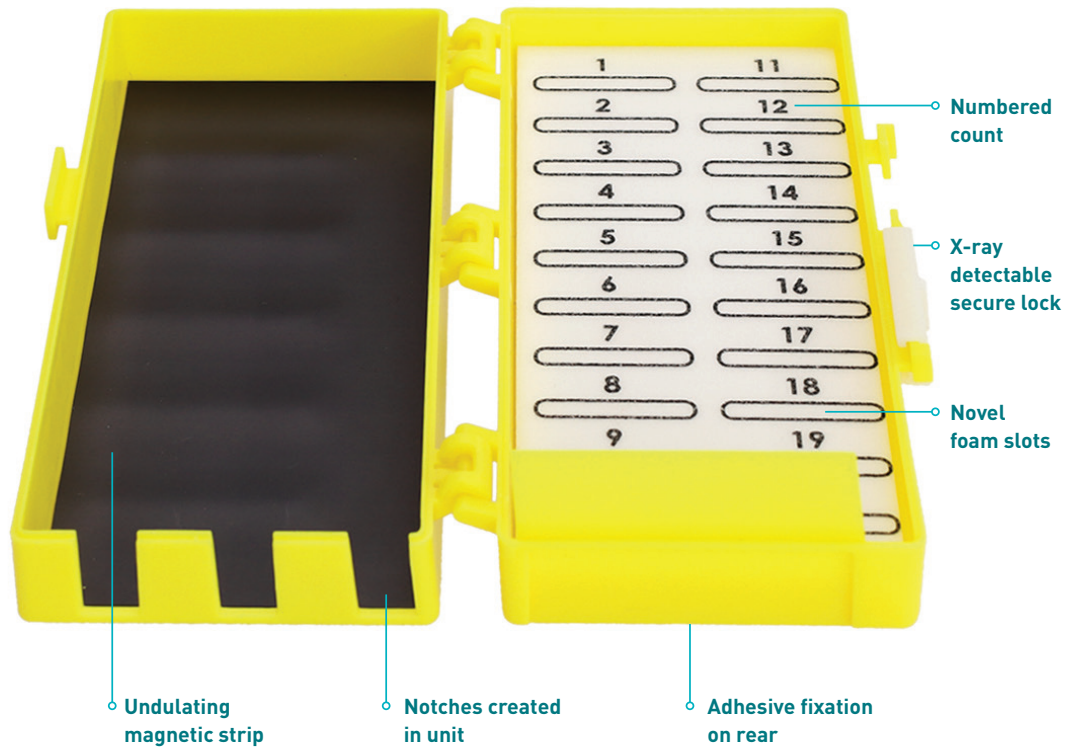

#### Sharps safety

Slots in the foam allow for easy insertion and removal of sutures without the need to reposition with hands or blunting the tip of the suture.



#### Easy blade lift

Undulating magnetic strip assists with blade retrieval from magnet.



#### Instrument organisation

Notches at one end of the unit can be used to rest instruments. The magnets or foam blocks can be used to secure sutures.

CLIMATE NEUTRAL<sup>1</sup>




As part of our commitment to better product stewardship, all SafeCount™ needle counters sold to Australian customers are supplied Climate Neutral.

## INSTRUCTIONS FOR USE

### HOW TO PARK AND REMOVE A SUTURE FROM THE FOAM SLITS



**01** Push the suture down into a foam slot and release from needle holder

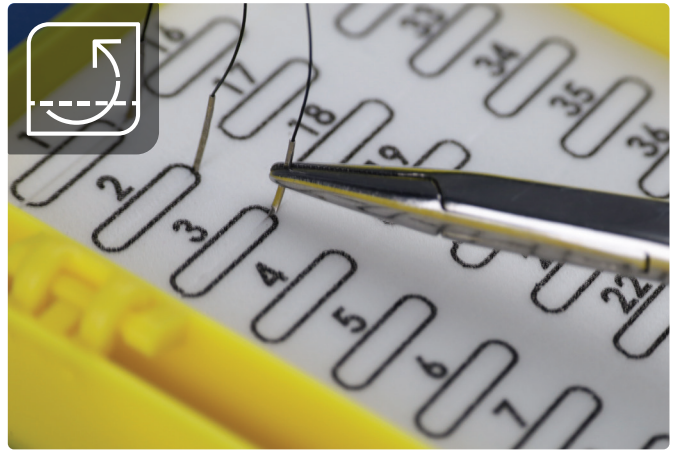

**02** Use the needle holder to pivot the suture so that the sharp end is fully covered below the surface of the foam



**03** To remove, secure base of suture with needle holder



**04** Extract suture

## PRODUCT DETAILS

<b>Description</b>	Needle counter 20 count/40 capacity, with foam block and magnetic in deep option
<b>Material</b>	Injection moulded polystyrene, magnet sheet, PE foam and paper
<b>Latex free</b>	Yes
<b>Shelf life</b>	—
<b>Non-Sterile code</b>	CEA-0010

<b>Size when closed</b>	37 x 127 x 74 mm
<b>Units per box</b>	20
<b>Boxes per carton</b>	4
<b>Carton quantity</b>	80
<b>Sterile</b>	Sterile
<b>Use</b>	Single use
<b>Colour</b>	Yellow, universal colour to indicate sharps

**Disclaimers:** 1. Multigate only purchases carbon credits via Carbon Neutral in approved, registered projects which have been independently assessed and verified as creating additional resources that help combat climate change using avoidance or sequestration means.

**IP Protected Product:** Intellectual Property rights such as patents, patent applications, registered designs, trademarks, etc. apply to this product group. See website for details.

**FOR MORE INFORMATION** [multigate.com.au](https://multigate.com.au) or call 1800 023 420 (Australia)  
[multigate.co.nz](https://multigate.co.nz) or call 0800 880 013 (New Zealand)