

Sterilisation Apparel & Accessories



High Performance Sterilisation Wrap
Innovation Focus: Barrier Performance

OCTOBER 2024

At Multigate we combine our deep understanding of how you work with our design capabilities to innovate and better serve you, your patients and our partners.

For 30 years we've created medical products, services and experiences that give you the advantage to perform at your best. We continue to do this.

Inspiring performances and better patient outcomes that will lead the world of hospital care to new standards.

We are Multigate.

We understand that sterile management, clinician protection and streamlined daily workflows are paramount in the CSSD.

With that in mind, we carefully designed our range of protective apparel, towels and tray liners and sterilisation accessories to power your performance.



MULTIGATE



Addressing CSSD clinician concerns about the time, costs and effectiveness of sterilising instruments, we introduced our simple to use, cost saving Pro-Zone® range of single use sterilisation wraps.

Made from polypropylene and using our 5 layer SSMMS protection technology, Pro-Zone® wraps offer an economic and time efficient sterilisation solution.



High Performance Sterilisation Wrap

Innovation Focus: Barrier Performance



TECHNOLOGY



Our Pro-Zone® range is made from polypropylene, utilizing state of the art 5 layer (Spunbond-Spunbond-Meltblown-Meltblown-Spunbond) technology.

DESIGN

Our entire Pro-Zone® range of sterilisation wraps are a versatile high quality wrapping system that delivers many performance advantages including:

- Strength and resistance to punctures, tears and abrasions*
- Anti-static coating reduces static “shocks” to clinicians
- Low linting
- Soft and easy to fold
- Excellent drapability and low memory
- Excellent microbial barrier
- Effective sterilant penetration

* when correct handling techniques are used

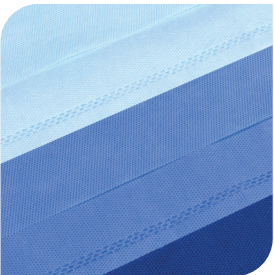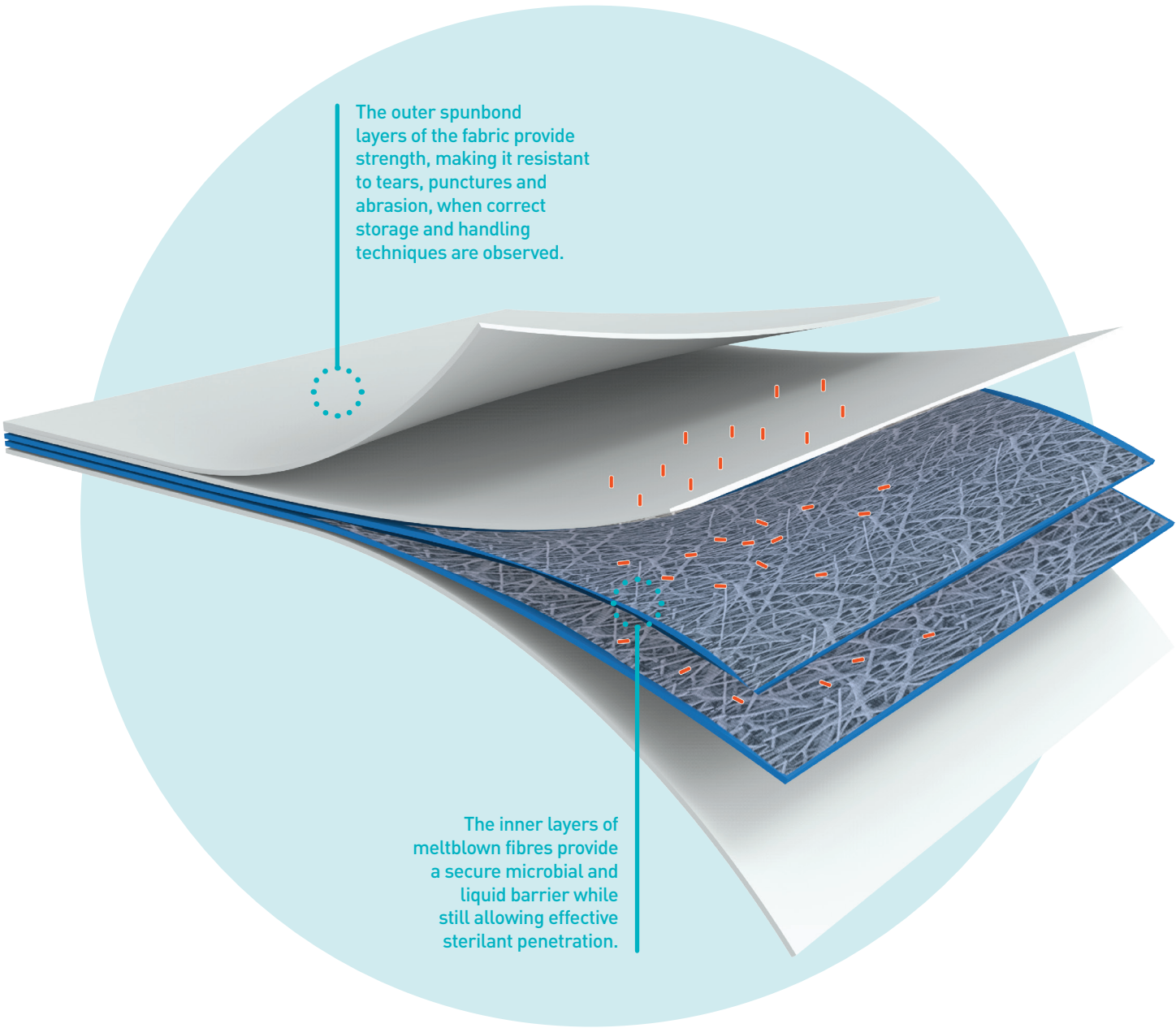

Pro-Zone®
The original Pro-Zone® sterilisation wrap.

Weights:
— 47gsm
— 61gsm
— 75 gsm

Compatible with:
— Steam
— ETO
— Low temperature gas plasma
— Low temperature steam/formaldehyde

Standards Compliance:
— AS/NZS 4187:2014
— ISO 11607



Pro-Zone®+2
Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.

Weights:
— 47gsm
— 61gsm
— 75 gsm

Compatible with:
— Steam
— ETO
— Low temperature gas plasma
— Low temperature steam/formaldehyde

Standards Compliance:
— AS/NZS 4187:2014
— ISO 11607



Pro-Zone®+4G
Pro-Zone®+4G offers the protection of 2 bonded sheets, plus the added advantage of an absorbent layer to wick any residual condensate that may compromise the integrity of the wrap.

Weights:
— 61gsm
— 75 gsm

Compatible with:
— Steam
— ETO
— Low temperature steam/formaldehyde

Standards Compliance:
— AS/NZS 4187:2014
— ISO 11607

INFECTION PREVENTION

Pro-Zone®’s strong 5 layer construction reduces the risk of tears, punctures and abrasions, and delivers a low linting, high performing sterile barrier system.

- ✓ Infection prevention
- ✓ Process efficiencies
- ✓ Comfortable handling
- ✓ Space saving
- ✓ Recyclable
- ✓ Low linting
- ✓ Anti-static

PERFORMANCE

Streamlined Process.

Our Pro-Zone® range eliminates the costly decontamination processes associated with reusable sterilisation systems and delivers faster instrument processing.

The fabric is soft, has low memory and provides excellent draping capability allowing Pro-Zone® wrap to double as a sterile field, permitting the aseptic removal of pack contents and further streamlining theatre processes.

In-Service Support.

We understand that every CSSD is different. Our Sterilisation Specialists are trained in CSSD principles and practises. They are available to:

- Provide advice
- Deliver in-service training and support
- Evaluate sterilisation wrap performance.

FEATURES

Comfortable Handling.

Pro-Zone® wraps are anti-static, easy to fold, gentle on your hands and light offering a packaging system that is easy to handle.

Space Saving.

Pro-Zone®’s compatibility with a broad range of sterilants and its ability to cater for all instrument shapes and sizes with just a few different weights and specifications, reduces the number of sterilisation components and returns valuable space to CSSDs.

Recyclable.

The entire Pro-Zone® range is recyclable.



Microbial barrier that delivers reliable barrier performance



Easy to fold material with low memory to provide excellent draping capability





Pro-Zone®

The original Pro-Zone® sterilisation wrap.



Pro-Zone®+2

Pro-Zone®+2 boasts the innovative use of 2 coloured sheets – blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap.

Pro-Zone® (47) Non-Sterile			COLOUR: BABY BLUE
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-465	300 x 300	47	1000
30-464	450 x 450	47	500
30-460	610 x 610	47	500
30-470	710 x 710	47	300
30-461	915 x 915	47	150
30-467	1000 x 1060	47	150
30-466	1000 x 1400	47	100
30-462	1215 x 1215	47	100
30-468	1300 x 1525	47	100
30-463	1830 x 1370	47	100

Pro-Zone® (61) Non-Sterile			COLOUR: BLUE
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-471	710 x710	61	200
30-472	915 x 915	61	150
30-475	1060 x 1200	61	50
30-476	1060 x 1400	61	50
30-473	1215 x 1215	61	50
30-479	1215 x 1215	61	50
30-477	1300 x 1525	61	50
30-478	1370 x 1370	61	50
30-474	1370x 1830	61	50

Pro-Zone® (75) Non-Sterile			COLOUR: DARK BLUE
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-484	1060 x 1200	75	50
30-483	1060 x 1400	75	50
30-480	1215 x 1215	75	50
30-481	1370 x 1830	75	40

Pro-Zone®+2 (47) Non-Sterile			COLOUR: BABY BLUE / LILAC
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-551Z	300 x 300	47	480
30-552Z	450 x 450	47	480
30-553Z	600 x 600	47	240
30-554Z	760 x 760	47	144
30-555Z	910 x 910	47	144
30-556Z	1010 x 1010	47	120
30-557Z	1210 x 1210	47	48
30-558Z	1370 x 1370	47	48
30-559Z	1370 x 1820	47	48

Pro-Zone®+2 (61) Non-Sterile			COLOUR: BLUE / LILAC
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-560Z	710 x710	61	96
30-561Z	910 x 910	61	72
30-562Z	1010 x 1010	61	72
30-563Z	1210 x 1210	61	24
30-564Z	1370 x 1370	61	24
30-565Z	1370 x 1820	61	24

Pro-Zone®+2 (75) Non-Sterile			COLOUR: DARK BLUE / LILAC
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)
30-570Z	910 x 910	75	72
30-571Z	1010 x 1010	75	48
30-572Z	1210 x 1210	75	24
30-574Z	1370 x 1370	75	24
30-573Z	1370 x 1820	75	24
30-575Z	1500 x 1500	75	24



Pro-Zone®+2, 4G

Pro-Zone®+2, 4G offers the protection of 2 bonded sheets, the innovative use of 2 coloured sheets, blue and lilac, bonded together, enabling clinicians to rapidly identify any tears or damage in the wrap. plus the **added advantage of an absorbent layer** to wick any residual condensate.

Pro-Zone®+2, 4G (61) Non-Sterile				COLOUR: BLUE
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)	
30-881	1010 x 1010 (absorbent layer 70 x 70)	61	24	
30-882	1210 x 1210 (absorbent layer 70 x 70)	61	24	
30-883	1370 x 1370 (absorbent layer 70 x 70)	61	24	

Pro-Zone®+2, 4G (75) Non-Sterile				COLOUR: DARK BLUE
RE-ORDER	SIZE (mm)	WEIGHT (gsm)	CARTON QTY (Pieces)	
30-891	1010 x 1010 (absorbent layer 70 x 70)	75	24	
30-892	1210 x 1210 (absorbent layer 70 x 70)	75	24	
30-893	1370 x 1370 (absorbent layer 70 x 70)	75	24	
30-894	1370 x 1820 (absorbent layer 70 x 70)	75	24	



Towels and Tray Liners

Made from non-woven spun lace fabric cellulose/ polyester (MediClean™) & cellulose/polyester/ viscose (MediCEL®). These towels are ideal for a wide range of dry or wet CSSD applications. They are ideal for use as a tray liner and applications where the towel is soaked in water for long periods of time.

MediClean Non-Sterile					
RE-ORDER	SIZE (cm)	COLOUR	PIECES/BOX	BOXES/CARTON	CARTON QTY (Pieces)
57-112	31.8 x 33 (small)	○ White	100	12	1200
57-113	60 x 35 (large)	○ White	100	12	1200

MediCel Non-Sterile					
RE-ORDER	SIZE (cm)	COLOUR	PIECES/BOX	BOXES/CARTON	CARTON QTY (Pieces)
57-120	35x30 (small)	○ White	100	12	1200
57-121	60 x 35 (large)	○ White	100	12	1200

Protective Apparel and Eyewear

- Thumb Secure Isolation/Protective Gown
- Liteback Thumb Secure Protective Gown
- Protective Eyewear

Thumb Secure Isolation/Protective Gown Non-Sterile					
RE-ORDER	COLOUR	SIZE	PIECES/BOXES	BOXES/CARTON	CARTON QTY (Pieces)
29-285	Yellow	Regular	15	5	75
29-286	Yellow	X-Large	15	5	75
29-385	Blue	Regular	15	5	75
29-386	Blue	X-Large	15	5	75

Liteback Thumb Secure Protective Gown Non-Sterile					
RE-ORDER	COLOUR	SIZE	PIECES/BOXES	BOXES/CARTON	CARTON QTY (Pieces)
29-291	Yellow	Regular	10	6	60
29-293	Yellow	XX-Large	10	6	60

Protective Eyewear Non-Sterile				
RE-ORDER	COLOUR	PIECES/BOXES	BOXES/CARTON	CARTON QTY (Pieces)
54-100NS	Black	15	10	150



Thumb Secure Protective Gown



Liteback Thumb Secure Protective Gown



Protective Eyewear



Protective Corner

Our protective corner helps to reduce piercing risk while wrapping.

Protective Corner Non-Sterile		
RE-ORDER	COLOUR	CARTON QTY (Pieces)
DP 2704	Orange	300

Empowering
performance
since 1986

Multigate Medical Products

27 Llewellyn Avenue

Villawood NSW 2163

Australia

T +61 (0)2 8718 2888

F +61 (0)2 8718 2889

enquiries@multigate.com.au

www.multigate.com.au

V.14.0

For more information please contact your
local Business Development Manager:

(02) 8718 2888 or visit **multigate.com.au**